# Ortho**P**revent**A**ligner<sup>®</sup>

OB-**M** soft

**D**

**GB**

**F**

**I**

**E**

## Patienteninformation

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D
<b>GB</b>
F
I
E

**Patient information**

## Prevention-oriented early treatment

Harmful factors that have a lasting effect on facial growth and tooth positions require **early preventive = precautionary** treatment. Harmful factors are sucking habits on the soother, thumb or other bodies as well as dysfunctions of the oral musculature, which are often associated with general postural weaknesses.

Harmful habits and dysfunctions such as the open mouth posture, incorrect tongue rest positions, muscle-weak lip closure as well as infantile swallowing (baby swallowing) are remedied with the **OrthoPrevent Aligner® OB-M**. In particular, frontal open bites, which result in speech disorders and restricted biting and chewing, can be corrected with the **OrthoPrevent Aligner® OB-M**. If this has resulted in a gaping tooth position with protruding teeth, these are encouraged to heal themselves. Regular wearing of the **OrthoPrevent Aligner® OB-M** is always a prerequisite.

## How jaw anomalies develop

Persistent **sucking habits** - of whatever kind - cause the tongue to be displaced backwards and downwards. It no longer rests against the teeth, thus causing an imbalance between the tongue and cheek muscles.

The upper jaw is constricted. There is a lack of space, which results in crooked incisors.

Even with an **open mouth posture**, the tongue only lies in the lower jaw and as such, it also causes the upper jaw to be too narrow.

When swallowing incorrectly, which also happens unconsciously about 1,500 times a day, the tongue always pushes itself between the incisors and shortens them.

The result: **An open bite.**

For the treatment of faulty swallowing and speech errors, logopaedic therapy should also be provided.

## The OrthoPreventAligner® OB-M

The **OrthoPreventAligner® OB-M** (OB = Open bite malocclusion) for the treatment of the frontal open bite is made of a highly elastic silicone. The lip shield for the upper jaw is extended and strongly tilted backwards. As a result, the shortened incisors and the curved jaw section are fully enclosed.

The effect of the **OrthoPreventAligner® OB-M** in eliminating the open bite is **indirect**: The tongue is prevented from permanently positioning itself between the incisors. This allows the teeth to lengthen again in a self-regulating manner. In addition, there is a **direct** effect by placing the upper lip shield on the usually protruding and shortened teeth, which are thereby directed towards the palate and simultaneously lengthened.

By closing the open bite, the defective phonation is corrected, pronunciation is improved and biting with the incisors is possible again.

**OB-M Deciduous dentition**

In the case of an open bite with prolonged vestibular pelotte  
*(Similar to Fig.)*


## About habits and malfunctions

The main causes for the development of open bite malocclusion are **sucking habits** on fingers, thumbs, soothers and other foreign bodies. The resulting jaw anomaly causes speech disorders (lisp) and restricts the bite function.

Open bite malocclusion is rarely inherited or congenital. It is an acquired (self-inflicted) "dynamic dysfunction" that can be corrected by the **OrthoPrevent Aligner® OB-M**. To make it easier to stop the harmful habits, the device also serves as a substitute for suckers or thumbs.

In addition, the **OrthoPreventAligner® OB-M** can successfully treat incorrect swallowing (a habit from infancy), the wrong tongue rest position, the open mouth posture and mouth breathing as "static functional disorders".

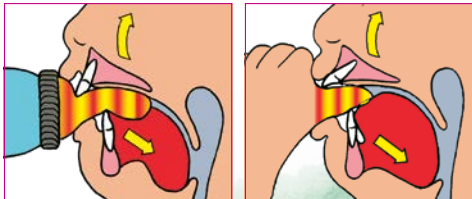By wearing the **OrthoPrevent Aligner® OB-M** regularly, these disorders and developed tooth misalignments are eliminated and normal oral functions are restored.

## Recommended wear

In order to break harmful habits and to eliminate **mal-functions** as well as to classify "crooked" teeth, the **OrthoPreventAligner®** should be worn regularly for at least 2 hours a day and at night.

## Wearing check

	Week 1		Week 2		Week 3	
	 hr.	 hr.	 hr.	 hr.	 hr.	 hr.
MON						
TUE						
WED						
THUR						
FRI						
SAT						
SUN						

## Instructions for use

Rinse the **OrthoPreventAligner®** with warm water before using for the first time. The **OrthoPreventAligner®** can be disinfected briefly by boiling water after infections.

## Cleaning instructions

Clean thoroughly with a toothbrush under warm running water. Alternatively, OROFAN® dental cleaning tablets can be used.

The **OrthoPreventAligner®** is not dishwasher or microwave safe. Only clean the hygienic **OrthoPrevent Aligner®** box with a damp cloth.

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